

Features

- Autonomously functioning cart with integrated silent oil-free compressor and spray water tanks
- High quality stainless steel housing
- 1 Three-way syringe (air/water/spray) with sterilisable tips
- Integrated spray water system for syringe and instruments (check valve prevents the back flow of water)
- 2 Spray water tanks: 2 × 1 litre
- 4 Instrument hoses with light (6-hole connector)
- Instrument hoses detachable
- Spray adjustment regulators for each instrument
- Foot control with spray on/off function and chip
- Integrated compressor: Capacity 75 l/min at 4 bar; Compressed air tank: 4 litres

Technical Data

- Air consumption instruments: 35-55 l/min
- Operating pressure handpieces: 2.5-3.5 bar
- Each instrument is individually adjustable
- Total weight: 62 kg
- Dimensions: 355 mm × 690 mm × 895 mm

Warranty: 2 years

Options: Electric motor with LED light, electric implantology motor with LED light, touchscreen and 99 programs per user, piezoelectric scaler, air-driven scaler, peristaltic pump, polymerisation lamp, turbine, air motor, etc.

P3-Cart

Freedom in its most elegant form


The P3-Cart combines maximum independence, uncompromising performance and an elegant, well-thought-out design – without any installation work. For commissioning, all you need to do is insert a single power plug – and the system is ready for use.

The handpieces are freely configurable and offer full performance even in continuous operation, meeting the highest professional standards. The cart remains pleasantly whisper-quiet and operates oil-free – for a calm, hygienic treatment environment.

The intuitive touchscreen allows every user to access up to 99 freely configurable programs. Depending on the configuration, treatments across all dental disciplines are possible – from general dentistry, prophylaxis, endodontics, orthodontics, and oral surgery to implantology.

The high level of operating convenience is also reflected in the externally removable handpiece hoses, practical connections for external devices and the dual spray water bottle solution: a simple switch allows seamless switching from one bottle to the other – without interrupting treatment.

The P3-Cart combines technical sophistication with elegant design – for dental professionals who expect flexibility, performance and top-level quality.

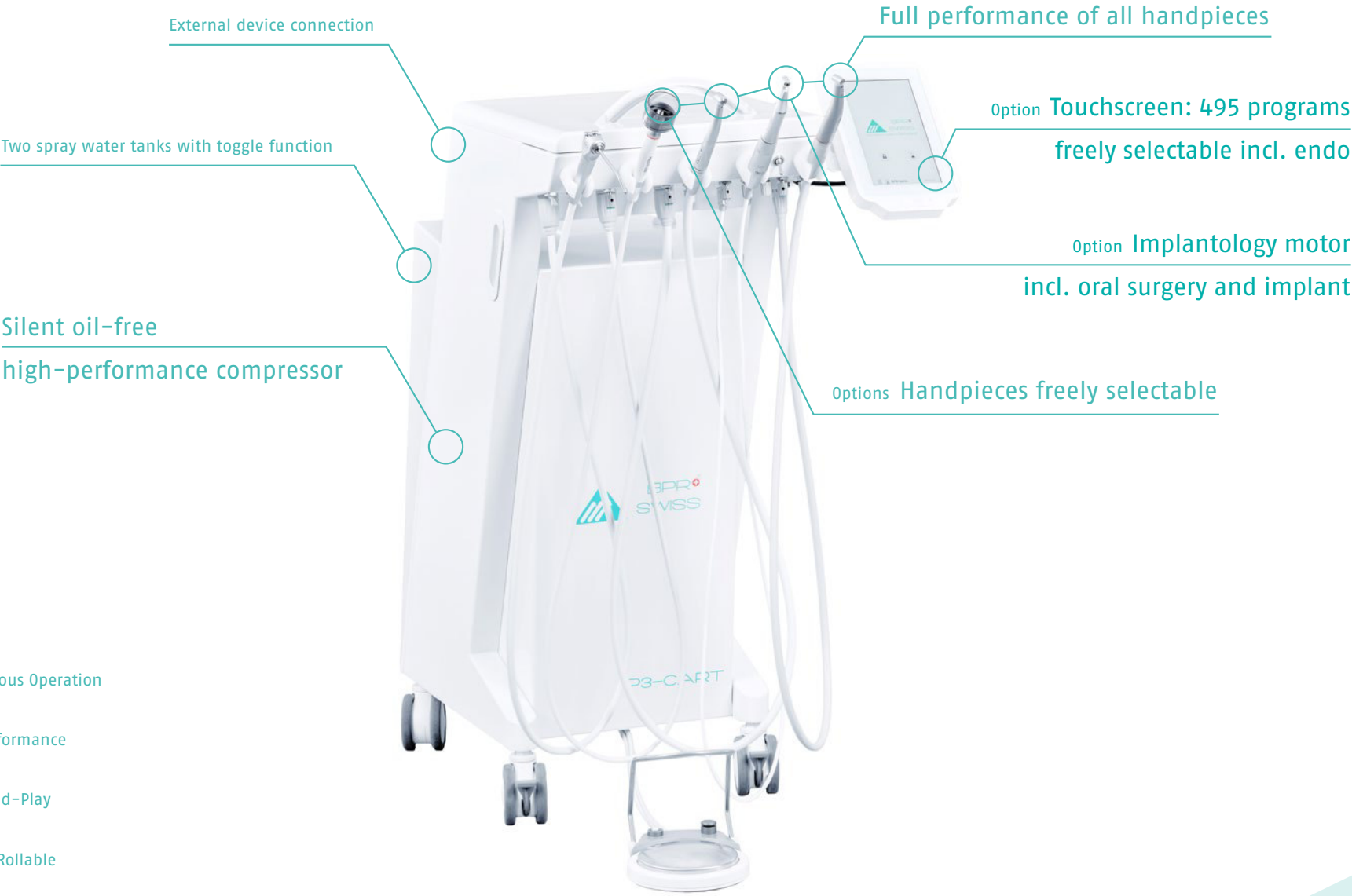

- Continuous Operation
- Full Performance
- Plug-and-Play
- Mobile/Rollable

P3-Cart

P3-Cart